



## Ultima Direct Puncture Drainage Set

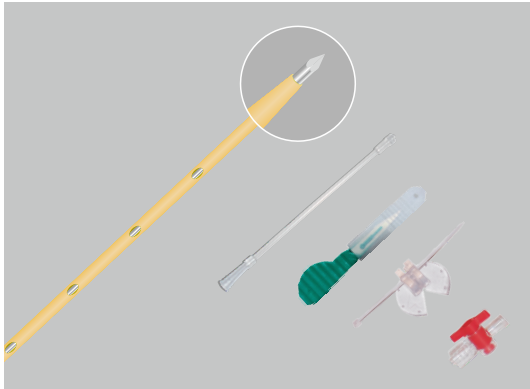

The Ultima Direct Puncture Drainage Set is used for percutaneous nephrostomy drainage following the Direct - Puncture technique, particularly when there is a dilated urine collection system or hydronephrosis.

- The catheter is back-loaded on a two-part trocar needle, straightened and then introduced directly through a skin incision
- The procedure does away with the use of guidewire and fascial dilators, as dilatation is not required
- The set may also be used for cyst and abscess drainage

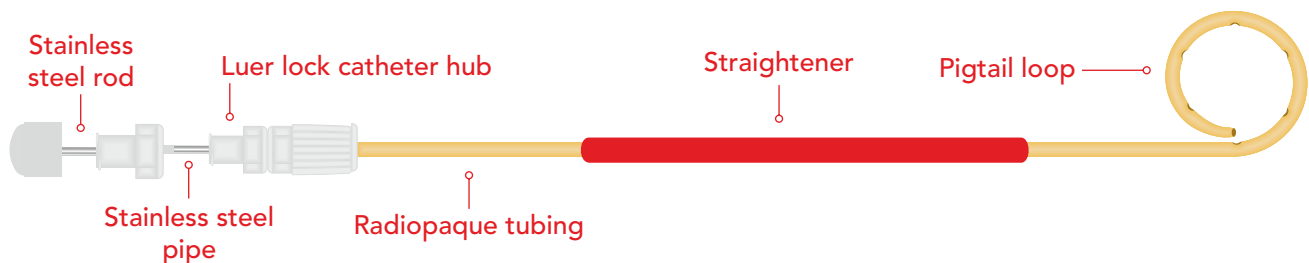

Order No.	Fr	Length (cm)	Qty. / Box
U-DPD-8/30	8	30	1
U-DPD-9/30	9	30	1
U-DPD-10/30	10	30	1
U-DPD-12/30	12	30	1
U-DPD-14/30	14	30	1

\*All products may not be available in all countries. Please contact our sales department for details.

### Set includes :

- Ultima Direct Puncture Drainage Set 30 cm
- Initial Puncture Needle Trocar Tip - 2 Part
- Urine Bag Connector
- Scalpel
- Catheter Fixing Disc
- Stopcock

### \*Combinations :

- With Locking Loop System : Add Suffix : LS
- With Hydrophilic Coating : Add Suffix : H