



## Trapper No Tip Nitinol Stone Basket

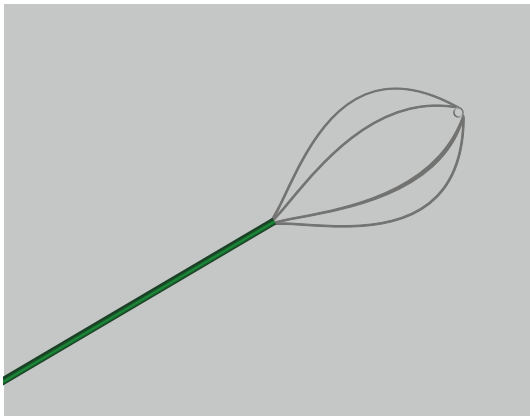

The Trapper No Tip Nitinol Stone Basket has been designed for endoscopically manipulating, dislodging, engaging and removing renal, ureteral or biliary stones.

- Tipless basket design to help with atraumatic retrieval of stones in the ureter and calyces
- Dismantleable handle
- Proximal end of our 1.8 Fr stone basket's shaft is made thinner for more flexibility
- Shaft and basket constructed of nitinol. The stone basket offers a kink-resistant and flexible design to help with stone extraction

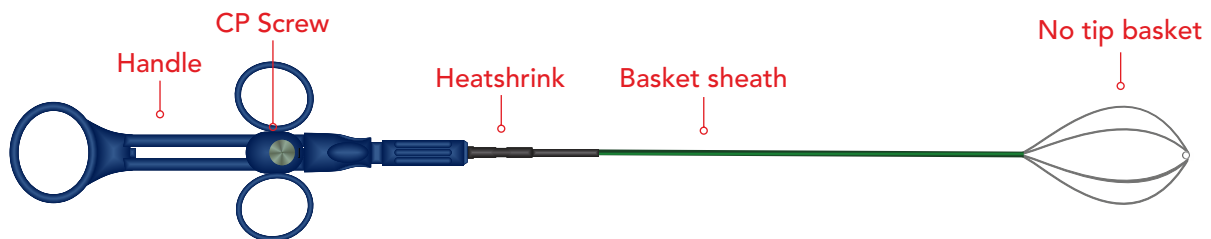

Order No.	Fr	Length (cm)	Qty. / Box
TNB4W/2.4F-90	2.4	90	1
TNB4W/2.4F-120	2.4	120	1
TNB4W/3F-90	3	90	1
TNB4W/3F-120	3	120	1
TNB4W/4F-90	4	90	1
TNB4W/4F-120	4	120	1

Order No. (with Polyimide Sheath)	Fr	Length (cm)	Qty. / Box
TNB4W/1.8F-90	1.8	90	1
TNB4W/1.8F-120	1.8	120	1

\*All products may not be available in all countries. Please contact our sales department for details.